


# LIMB-O™ ... When Less Is More

The innovative streamlined design of the Limb-O makes it the ideal breathing circuit for ventilation.

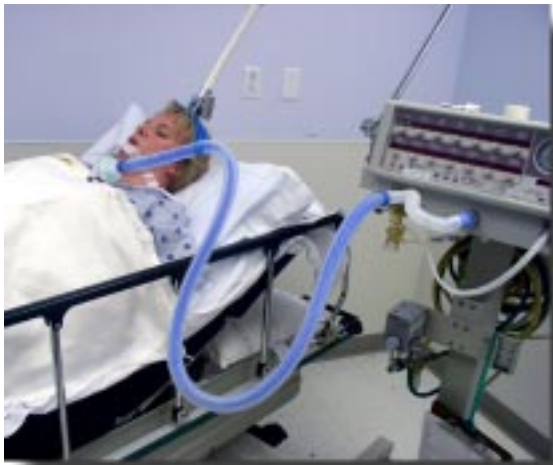

**LESS** Weight = *Limb-O* weighs 25% less than traditional two-limb circuits which reduces torque on the patient's airway.

**LESS** Inventory = *Limb-O* consolidates your inventory! The size of the lumen accommodates both adult and pediatric patients.

**LESS** Compliant = *Limb-O's* low compliance tubing allows for delivery of a more accurate tidal volume.

**LESS** Heat Loss = *Limb-O's* common wall/septum separating inspiratory and expiratory lumens permits heat transfer from warm exhaled gasses.



*Limb-O* offers **MORE** versatility with this single limb circuit when used as a transport circuit and also in the OR.



20 Campus Road, Totowa, NJ 07512, USA • TEL: 800-932-0760 • FAX: +1-973-790-3307 • www.vital-signs.com  
Sussex Business Village, Lake Lane, Barnham, PO22 0AL, UK • TEL: +44-1243-555300 • FAX: +44-1243-555400 • vitalcare@vital-signs.co.uk

vital signs<sup>in</sup>

vital

# LIMB- $\Theta$ ™



<b>Product Description</b>	<b>Limb-<math>\Theta</math>™ Adult and Pediatric Volume Ventilator Circuit</b>												
<b>Ordering Information (20/case)</b>	<b>Product Code</b>	<b>Description (Code assignment depends on configuration; basic configurations are as follows)</b>											
	AFNXXXX X	72" length											
<b>Accessories</b>	5510	Limb- $\Theta$ Transport System - 7' Star Lumen Oxygen Tubing set and APL Valve Set											
<b>Size</b>	72"												
<b>Weight with Connectors<sup>1</sup></b>	170g												
<b>Compliance of Circuit @ 100 cm H<sub>2</sub>O<sup>1</sup></b>	.93 ml/cmH <sub>2</sub> O												
<b>Insp. &amp; Exp. Resistance at Various Flow Rates<sup>1</sup> (cmH<sub>2</sub>O)</b>	Length	10 L/min		20 L/min		30 L/min		40 L/min		50 L/min		60 L/min	
		I	E	I	E	I	E	I	E	I	E	I	E
	72"	.14	.16	.46	.57	1.06	1.24	1.80	2.12	2.74	3.28	3.69	4.39
<b>Circuit Materials</b>	Corrugated Tubing and Septal Wall - Polyethylene/EVA 22mm Connectors – Medical Grade Polypropylene Cuff Band - Polyisoprene Expandable Extension Hose - Polypropylene												

20 Campus Road, Totowa, NJ 07512, USA • TEL: 800-932-0760 • FAX: +1-973-790-3307 • www.vital-signs.com  
 Sussex Business Village, Lake Lane, Barmham, PO22 0AL, UK • TEL: +44-1243-555300 • FAX: +44-1243-555400 • vitalcare@vital-signs.co.uk

**vital signs inc.**

<sup>1</sup>Test data on file.